




Although our packaging is designed to be strong and durable, please handle and unpack the package carefully to avoid accidental damage to its contents.

Need support? Feel free to reach out to us at customer.service@quantoom.com

Instructions for Safely Unpacking the Ntensify[®] micro

STEP 1: Opening of the wooden box

You will receive the **Ntensify[®] micro** in a fully closed and secured wooden box designed to protect the device as effectively as possible against potential impacts during transport. To open the box, use a hammer or a screwdriver to lift the two metal tabs located on the front side of the box. This step should be carried out by a specialist using the appropriate tools.

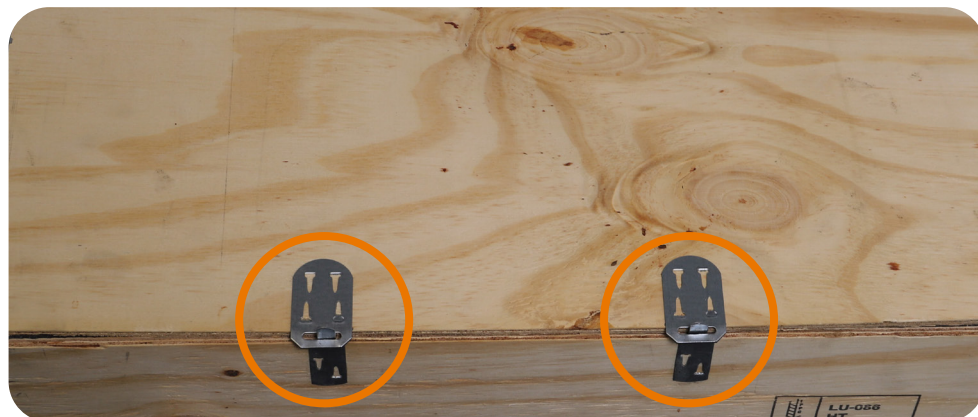


STEP 2: User Manual and Quick Start Guides

After opening the box, you will find several official documents related to the **Ntensify[®] micro**.

- On the left:
 - The User Manual
- On the right:
 - The present Unpacking process guide
 - The Quick Start Guide – 20 mL configuration
 - The Quick Start Guide – 96-well configuration
 - A marketing Guide


Please make sure to read these documents carefully before using the **Ntensify[®] micro**.

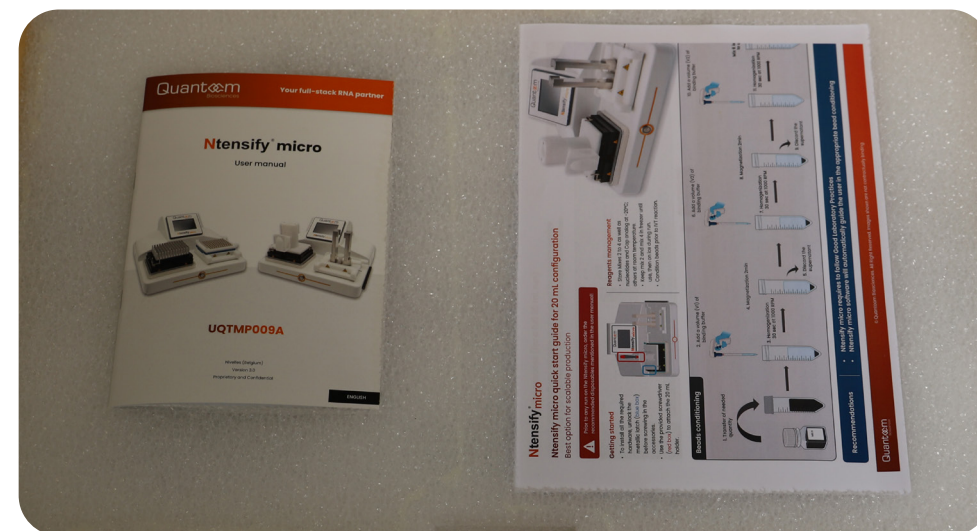

STEP 3: Accessories Overview

Once the first foam layer has been removed, you will find the accessories required to operate the **Ntensify[®] micro**:

- A power cable (n° 1)
- A touch stylus (n° 2)
- Angular adaptor (n° 3)
- 20 mL agitation accessory (n° 4)
- Screwdriver (n° 5)
- 20 mL magnetization accessory (n° 6)
- 96-well agitation accessory (n° 7)
- 96-well magnetization accessory (n° 8)

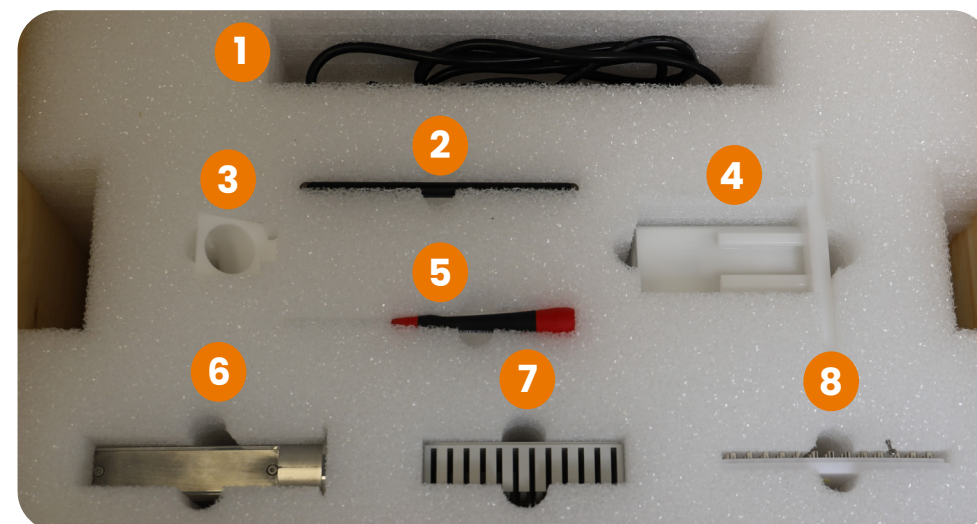


STEP 4: Protection for the Ntensify® micro

Once the second foam layer has been removed, you will see the **Ntensify® micro**, protected by a third foam layer and a lid stabilizer. You can carefully remove this third layer.



STEP 5: Removing of the stabilizer

Simply remove the yellow stabilizer used to protect the lid of the **Ntensify® micro**.

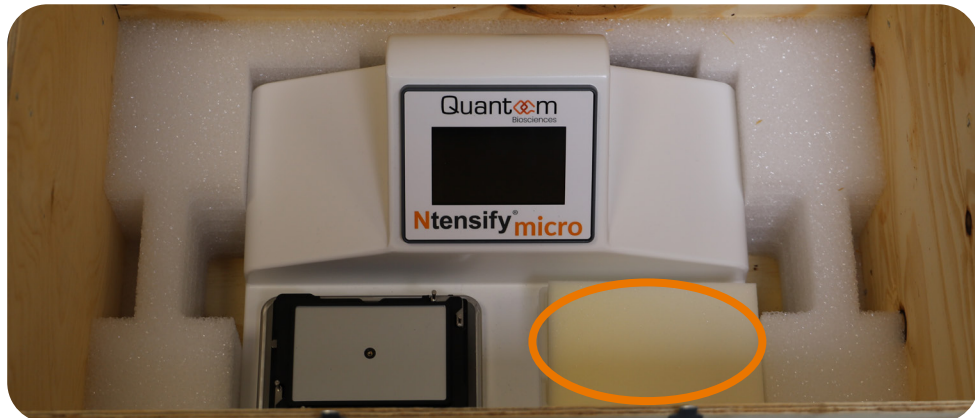

STEP 6: Removing of the Ntensify® micro

You can now safely handle the **Ntensify® micro** either by using the internal handles directly attached to the device or by using the external handles to remove the fourth foam layer, which will further facilitate unpacking. Be careful: the device weighs around 15 kg (33 lbs).



STEP 7: Box and foam layers management

After removing the **Ntensify® micro**, please place all foam layers back into the box and store it safely in case the device needs to be transported in the future.

