
Set-up & beads conditioning

Ntensify® micro

⚠ When using the Ntensify® micro, order the following disposable: 50 mL centrifuge tube DNase/RNase free.

Getting started

1. Unlock the metallic latch (**in blue**) **before** screwing in the accessories.
2. Attach the 20 mL holder with the provided screwdriver (**in red**).

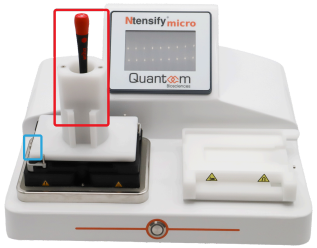


Run preparation

- Thaw all mixes at room temperature.
- Keep Mix 2 and Mix 4 in the freezer until use.
- Prepare the linearized DNA at 0.25 mg/mL for saRNA.
- Prepare the necessary volume of Ethanol 80%.

Beads conditioning

- Must be performed before each production.
- Conditioned beads must be used within 24h after preparation.

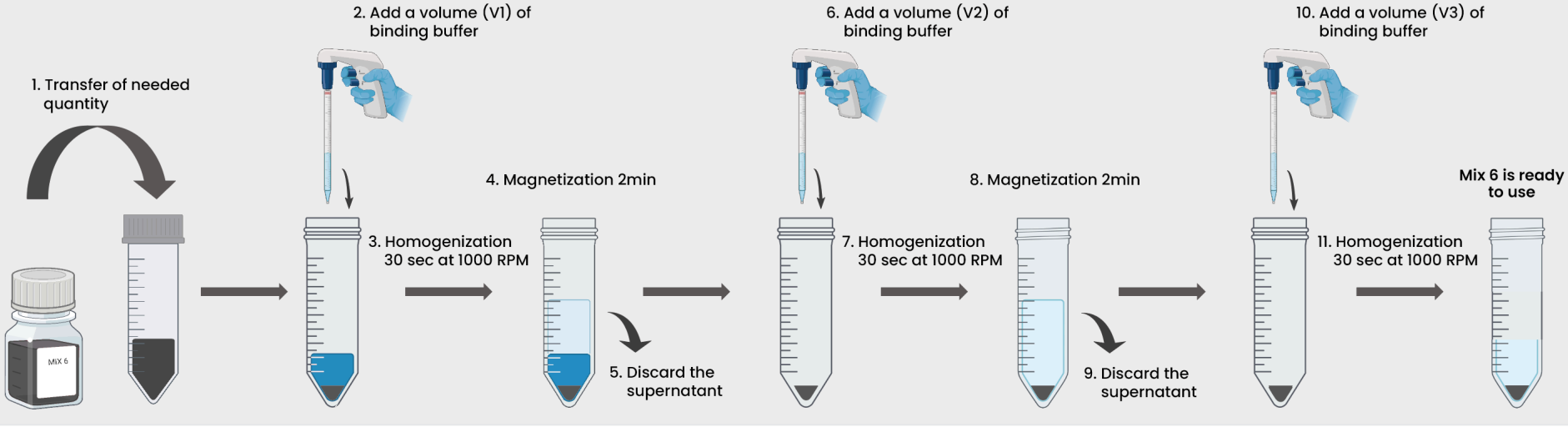

RECOMMENDATIONS

Follow Good Laboratory Practices at every step of the run.
The software guides you through all the production steps.

© 2026 Phoenix Biosciences SA, d/b/a Quantoom Biosciences. All rights reserved. Images shown are not contractually binding.


Overview of the process

Ntensify® micro

The buzzer rings after incubation, centrifugation and at the end of the run; respect the timing given by the software.

USP

I. Pre-mix heating

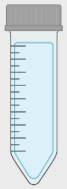


Add in the following order:

1. Mix 3
2. Selected uridines (rATP, rUTP, rGTP, rCTP)
3. Cap Analog
4. Add mix 2 to the previous premix volume.

Incubation · min. 15 min


II. IVT

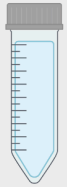

5. Add linearized DNA

Min. 60 min


III. DNA removal

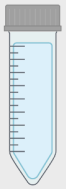

6. Add mix4

Min. 10 min

IV. Quenching

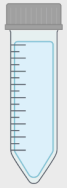

7. Add mix 5


Min. 15 min

DSP


REPEAT 2 TIMES

V. Binding


8. Add conditioned beads

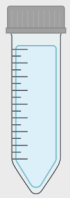

9. Add Ethanol 100%

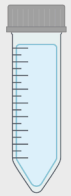

Magnetic separation

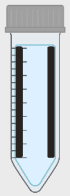


Removal of supernatant

VI. Washing

10. Add Ethanol 80%



Magnetic separation

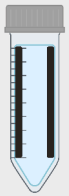


Removal of supernatant


Wash 4 times

VII. Elution

11. Add RNAse free water

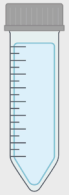

Magnetic separation

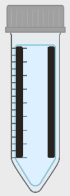

Collect the purified RNA

Elute 2 times

RECOMMENDATIONS

During magnetization, the tube can be tilted to facilitate supernatant pipetting.

The average duration of a run on the Ntensify® micro in the 20 mL configuration is about 4 hours (including preparation).

© 2026 Phoenix Biosciences SA, d/b/a Quantoom Biosciences. All rights reserved. Images shown are not contractually binding.